

LIFE BP Pro

Wireless Blood Pressure Monitoring System

The LifeStar™ BP Pro wireless blood pressure and pulse rate monitor is an advanced device for measuring and assessing hypertension. Monitoring of blood pressure can help avoid long term complications of cardiovascular diseases.

The LifeStar™ BP Pro communicates with the LifeStar™ handheld wireless device and the LifeStar™ Network, provided by LifeWatch Technologies.

- Wireless non-invasive blood pressure and pulse rate monitor
- Portable and simple to operate
- Comfortable and rapid testing
- Automatic inflation and deflation
- Anywhere, Anytime



Item #: FG-00017C




Item #: 903-7006-00C

After performing a test, the results are automatically transmitted via Bluetooth technology to the LifeStar™ handheld device, where it can be viewed and saved. When convenient, the data can be transmitted to the LifeStar™ Network, where physicians may view patient data. Additional remote patient monitoring services are also available. The LifeStar™ Network provides privacy and data protection of patient information, and has the ability to track vital signs data over time.



LIFE BP Pro

Wireless Blood Pressure Monitoring System

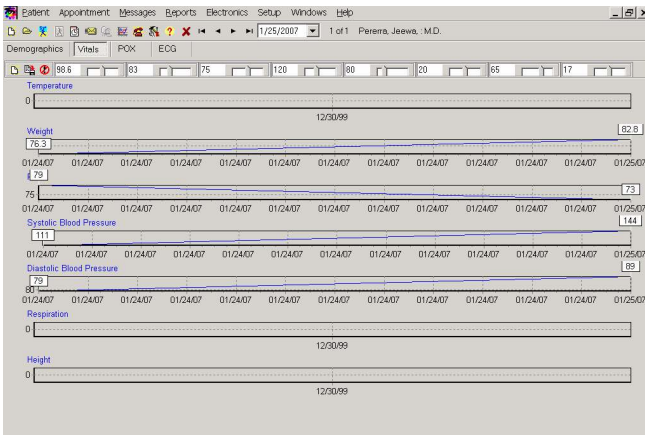

Technical Specifications	Parameter	Value
	Measuring principle:	Oscillometric
	Measurement range	
	Systolic (SYS)	70 - 260 mmHg
	Diastolic (DIA)	45 - 180 mmHg
	Pulse (Pulse)	40 - 240 1/min
	Display range:	0 - 300 mmHg/Pulse: 40 - 240 beats per minute (bpm)
	Storage capacity	99 measurements (systolic and diastolic pressure, pulse rate, date and time)
	Accuracy:	± 3 mmHg in display range
	Data security	Non-volatile memory
	Data transfer	Bluetooth technology
	Power source	2 alkaline 1.5V batteries (AA, Lr06)
	Measuring time	Approximately 30 seconds
	Controls:	2 multifunction keys: Start/Stop and Menu
	Display	LCD visible size approx. 45 x 40 mm
	Operating temperature	+10 °C to +40 °C
	Operating humidity	15% - 90%, non-condensing
	Storage temperature	-20 °C to +50 °C
	Storage humidity	15% to 95%
	Dimensions	125 mm W x 131 mm D x 51 mm H
	Weight	270 g



Distributed by LifeWatch Technologies Inc. ©2007. All rights reserved. LifeStar, the LifeStar symbol and HealthKit are trademarks or registered trademarks of LifeWatch Corp.

LifeWatch Technologies Inc.
 O'Hare International Center II
 10255 West Higgins Road Suite 120
 Rosemont, IL 60018
 Toll Free: 800-633-3361
 Fax: 800-954-2375
www.lifewatchtechnologies.com