

The Simplest Method For Monitoring Your Pacemaker Patients In The Convenience Of Their Home

The Pacemaker Monitoring System *CarryAll® Monitor*

The CarryAll® pacemaker monitor provides a simple, effective solution for monitoring your pacemaker patients in the convenience of their home.


The portable, battery-operated CarryAll monitor transmits real-time ECG/pacemaker activity along with measured atrial and ventricular pacemaker pulse widths over the telephone.

Digital signal processing achieves very high performance in discriminating between pacing pulses and non-pacing signals (including transthoracic bio-impedance pulses.)

Color-and shape-coded switches, large text instructions and pictures encourage easy step-by-step compliance by patients. The Patient Alert function provides both visual and audible indicators. The battery can be easily replaced and an automatic shut-off feature prevents unnecessary battery depletion.

The CarryAll monitor includes wrist electrodes, a magnet and a built-in telephone cradle modem in a single, easy-to-carry unit.

- Compatible with Single Chamber, Dual Chamber and Rate-Responsive Pacemakers in Unipolar or Bipolar Pacing Modes
- Digital Design Provides Accuracy and Reliability
- Selectable Transmission Modes Controlled by the Clinician
- Easy-to-Use Design Enhances Patient Compliance

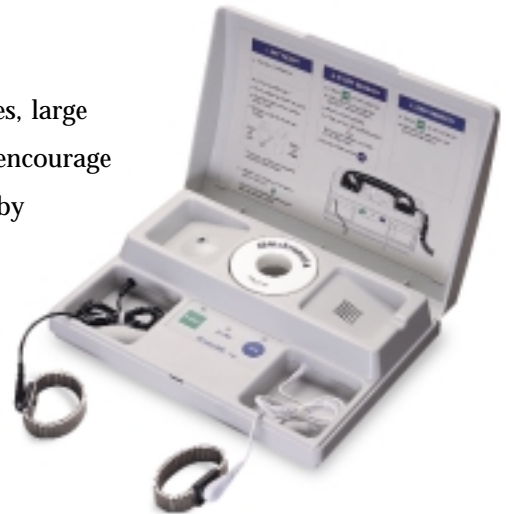

Instromedix
A Card Guard® Company

The Pacemaker Monitoring System

CarryAll[®] Monitor



Transmission Specifications	Models	PDT	DDD	
	Pulse Width Data Mode	■	■	
	ECG Mode	■		
	Enhance Mode	■		
	Center Frequency (typical)	1965 Hz	1965 Hz	
	Deviation	100 Hz/mV	100 Hz/mV	
Physical	Dimensions	7.5 x 10.5 x 2.3 in 19.1 x 26.7 x 5.8 cm	Weight with magnet and battery	2.6 lb 1.18 Kgm
Electrical	Bandwidth	0.5 to 100 Hz		
	Input Impedance	2.4 Mohm		
	Pulse Width Range	0.2 to 2.50 msec		
	Pulse Width Accuracy	±0.02 msec		
	Pulse Sensitivity	1mV, 1 µsec rise and fall time		
	Refractory Period	40 msec		
Battery	Battery Type	9V Alkaline		
	Battery Life	10 hours (continuous operation)		

Instromedix
A Card Guard[®] Company

O'Hare International Center II
10255 West Higgins Road, Suite 120
Rosemont, IL 60018

t 800.633.3361
f 800.954.2375

www.instromedix.com

Copyright ©2003 Card Guard[®] Group of Companies. All rights reserved. Card Guard[®] reserves the right to change specifications at any time without notice. Card Guard, Instromedix, Instromedix logo, and CarryAll are trademarks or registered trademarks of Card Guard[®] Group of Companies. All other trademarks belong to their respective owners. No affiliation with these trademarks or their owners is expressed or implied herein.
IMX038 Rev A M07-0023